


# Operation Manual Model: 9660

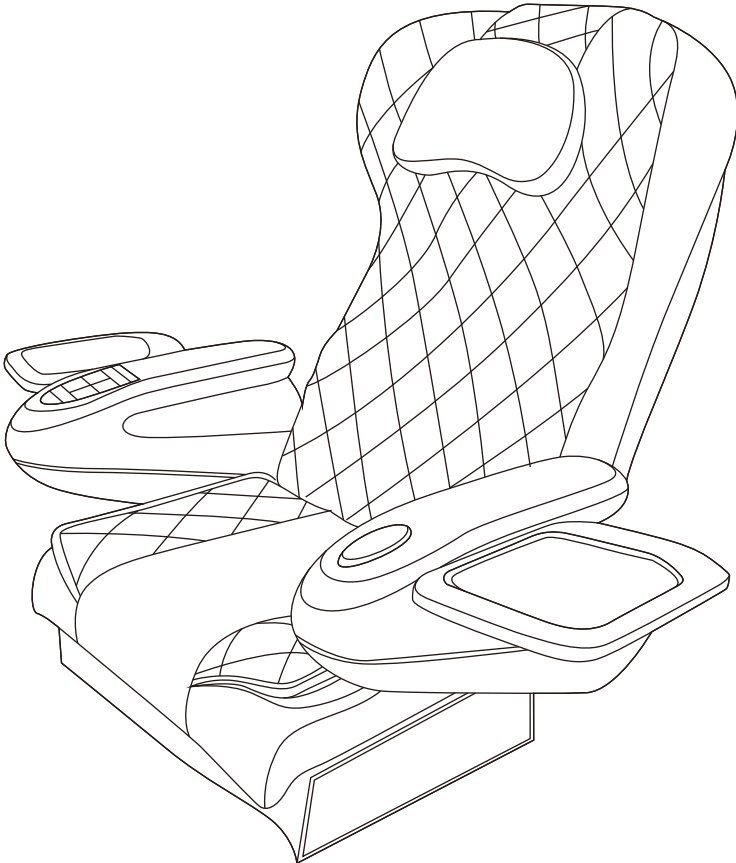

Please read this manual before operating

**9660 Massage Chair®**, **Gulfstream Inc.®** are trademarks of **Gulfstream Inc.**

Patent #D709619, #8726430, #8680699, #D701318, #8662848, #8531048, #8366418, #D670938, #D670354, #D662340, #D660444, #D650084, #D646082, #D636090, #D632798, #D628704, #D615326, #D613099, #D612509, #D610697, #D600950, #D596883, #D596882, #D590188, #D589285, #D588844, #D583960, #D583576, #D583171, #D565739, #D556480, #7111334, #6991182

© 2021 All rights reserved. Design and Specifications are subject to change without notice.

## **Welcome to True Touch® Massage Chair**

True Touch® chairs are the perfect bridge between state-of-the-art technology and the ancient, healing art of massage. True Touch® massage chairs utilize trademarked Gulfstream Inc. massage technology. This patented shiatsu system emulates the back and spinal care techniques used by massage and chiropractic professionals.

To ensure that you use this chair safely and properly; please read this manual, in its entirety, before beginning usage.

### **About Gulfstream Inc.**

Based in North America, Gulfstream Inc. is a leader in salon/spa equipment and furniture. Gulfstream's success is attributed to innovation and cutting edge products. Superior craftsmanship and ingenious design has led to some of the industry's most luxurious products. From custom color glass bowls for pedicure spas, to environmentally friendly, ultra hygienic poly liners for their pedicure bowls, Gulfstream Inc. leads the way for the salon/spa industry.

On-site production enables Gulfstream Inc. to continually upgrade their products, thus ensuring that their customers receive the latest designs and technology. Gulfstream Inc. is committed to meeting the evolving needs of the worldwide salon/spa market. Gulfstream sets the trends that others follow.

For more details and information about Gulfstream Inc. please visit our website at [www.gulfstreaminc.com](http://www.gulfstreaminc.com) or call our service centre at 1-877-622-4853.

Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our Company reserves the right to revise the design and description of this product without any further notice. Please refer to the actual color of this product.

## CONTENTS

---

Important Safety Instructions	03
Safety and Maintenance	04
Massage System Installation	05
Name and Function of Components	06
Function	07
Method of Usage	08-15
Product Specification	16

**Warning:** ⚠

This product can expose you to chemicals including lead which are known to the State of California to cause cancer or birth defects or other reproductive harms. For more information go to [www. P65Warnings.ca.gov/furniture](http://www.P65Warnings.ca.gov/furniture).

**When using an electrical appliance, basic precautions should always be followed, including the following:**

Read all instructions before using (this appliance)

**DANGER-**To reduce the risk of electric shock:

1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

**WARNING-**To reduce the risk of burns, fire, electric shock, or injury to Persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and Cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- 4) Use the appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) The appliance is intended for household and indoor use only.
- 14) Connect this appliance to a properly grounded outlet only ,see grounding instructions .
- 15) Keep children away from extended foot support (or other similar parts )
- 16) Do not use massager in close proximity to loose clothing or jewelry.
- 17) Keep long hair away from massager while in use.
- 18) Connect this appliance to properly grounded outlet only. See Grounding Instructions.
- 19) Donot use massager inclose proximity to loose clothing or jewelry.
- 20) Keep long hair away from massager while in use.

# SAVE THESE INSTRUCTIONS


## GROUNDING INSTRUCTIONS

This product must be grounded ,if it should malfunction or breakdown ,grounding provides a path of least resistance for electric current to reduce the risk of electric shock ,this product is equipped with a cord having an equipment grounding conductor and a grounding plug ,the plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.


**DANGER** :Improper connection of the equipment grounding conductor can result in a risk of electric shock ,check with a qualified electrician or serviceman of you are in doubt as to whether the product is properly grounded ,do not modify the plug provided with the product-if it will not fit the outlet,have a proper outlet installed by a qualified electrician.

1)This product is for use on nominal 120-volt circuit,and has a grounding plug that looks like the plug illustrated in sketch on the page 16 .A temporary adapter that looks like the adapter illustrated in sketches B and C may be used to connect this plug to a 2-pole receptacle as shown in sketch B if a properly grounded outlet is not available the temporary adapter should be used only until a properly grounded outlet (shetch A )can be installed by a qualified electrician ,the green colored rigid ear ,lug ,or the like extending from the adapter must be connected to a permanent ground such as a properly grounded outlet box cover ,Whenever the adapter is used ,it must be held in place by a metal screw .

Installation, maintenance, cleaning by professional personnel.

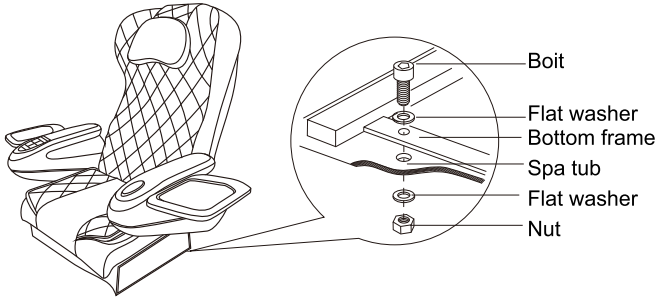

1. select a concrete or wood platform with an area more than 40cm\*36cm
2. drill four holes on each corner with a diameter 7.8mm
3. locked the massage chair on the platforms

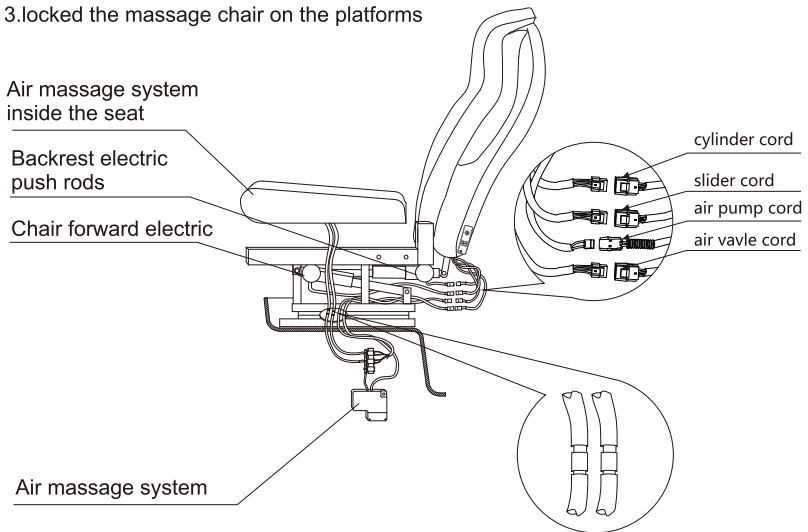


Fix it one by one the No

NOTE: Remark: air pump box(include air pump, magnetism valve and so on) is an optional parts. Clients can choose it at their request.

Connection of air pump.

1. Air pump box placed on a horizontal surface.
2. Connect the air pump cord to massage chair H to H, I to I.
3. Connect the magnetism cord to massage chair J to J, K to K.
4. Connect the slider and cylinder cord to massage chair M(M1) to M(M1), N(N1) to N(N1).
5. Connect black pipe to black pipe, connect white pipe to white pipe. Connect S1, S2 to massage chair relative pipes.

## ■ Name and Function of Components

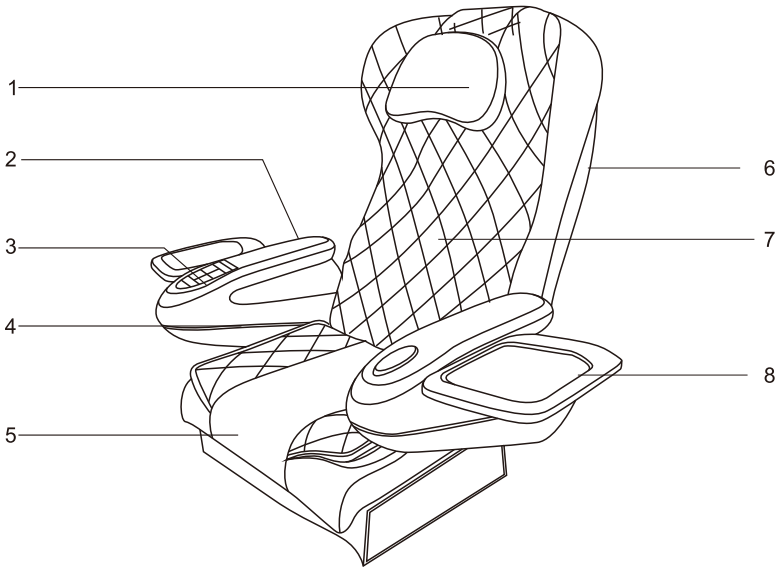

1. Detachable Pillow Pad
2. USB Charging Port
3. Controller
4. Armrest
5. Detachable seat cushion
6. The Back Covering
7. Back Pad
8. Tray

- Designed with a four-wheel driven muted massage mechanism which can move up and down in the backrest.
- Air pressure function: Buttock is equipped with 5 airbags, with 2 modes and 3 strengths.
- Programmed with five kinds of emulational massage functions of Shiatsu, kneading, flapping, knocking and synchronous kneading and flapping.
- Programmed with manual and automatic control modes for the upper body.
- Under the automatic mode, it offers three automatic changed combinations of massage function: lower, upper and full. Under the manual mode, it offers three options of overall, fixed and partial area. Under fixed state, the massage hands can be moved up and down; Each massage function has five levels of speed options; the width between the two massage hands can be adjustable in three levels of wide, middle and narrow when it is under flapping, shiatsu or knocking state.
- Seat frame forward and backward movement, backrest reclining and declining movement.
- Restoring the massage hands after turning off the machine.
- Design with USB charger on right armrest.



# Method of Usage



## 1. Controller operate

Button	Adjust massage balls
	To start or stop massage of the upper body

### (1) Auto Massage

Button	Adjust massage balls
	To carrying on automatic massage around the neck shoulder mainly
	Carry on automatic massage around in carrying the waist mainly
	Adopt the various tactics to carry on automatic massage to different position of the upper part of the body, the massage position covers the whole upper part of the body.

### (2) Upper body manual massage

Button	Function	Display status	Adjust massage balls
	Knead	<ul style="list-style-type: none"> <li>● Knead</li> <li>○ Flap</li> <li>○ Knock</li> <li>○ Press</li> </ul>	Speed
	Flap	<ul style="list-style-type: none"> <li>○ Knead</li> <li>● Flap</li> <li>○ Knock</li> <li>○ Press</li> </ul>	Speed, Width
	Knock	<ul style="list-style-type: none"> <li>○ Knead</li> <li>○ Flap</li> <li>● Knock</li> <li>○ Press</li> </ul>	Speed, Width
	Press	<ul style="list-style-type: none"> <li>○ Knead</li> <li>○ Flap</li> <li>○ Knock</li> <li>● Press</li> </ul>	Speed, Width
	Kneading & Flapping	<ul style="list-style-type: none"> <li>○ Knead</li> <li>○ Flap</li> <li>○ Knock</li> <li>○ Press</li> </ul>	Speed
	Message modes are equipped with 5 strengths.		

## Method of Usage



Intensity		For speed and massage treatments selection. five intensity options available. the indicators are as follows(black stands for on)	<table style="display: inline-table; vertical-align: middle;"> <tr> <td>● I</td><td>● I</td><td>● I</td><td>● I</td><td>● I</td></tr> <tr> <td>● II</td><td>● II</td><td>● II</td><td>● II</td><td>● II</td></tr> <tr> <td>● III</td><td>● III</td><td>● III</td><td>● III</td><td>● III</td></tr> <tr> <td>1</td><td>2</td><td>3</td><td>4</td><td>5</td></tr> </table>	● I	● I	● I	● I	● I	● II	● II	● II	● II	● II	● III	● III	● III	● III	● III	1	2	3	4	5
● I	● I	● I	● I	● I																			
● II	● II	● II	● II	● II																			
● III	● III	● III	● III	● III																			
1	2	3	4	5																			
Width		Width:when massage treatment of flapping,shiatsu,or knocking has been chosen,press this button to adjust the width of the massage hands,three options,not available under auto kneadingstate.	<table style="display: inline-table; vertical-align: middle;"> <tr> <td>● I</td><td>● I</td><td>● I</td></tr> <tr> <td>● II</td><td>● II</td><td>● II</td></tr> <tr> <td>● III</td><td>● III</td><td>● III</td></tr> <tr> <td>narrow</td><td>mid</td><td>wide</td></tr> </table>	● I	● I	● I	● II	● II	● II	● III	● III	● III	narrow	mid	wide								
● I	● I	● I																					
● II	● II	● II																					
● III	● III	● III																					
narrow	mid	wide																					
Back Stretching		Upper body(overall) trip massage																					
		Small area(partial) trip massage																					
		Fixed position (fixed) massage																					
Roller Position		When fixed massage has been chosen, press this button to move the massage device downward to suit precise massage points, not available under auto state.																					
		When fixed massage has been chosen, press this button to move the massage device upward to suit precise massage points, not available under auto state.																					





### (3)Other manual function

Timing		Three modes of timing: 10 mins, 20 mins and 30 mins.
--------	--	--

### (4)Air mode


Air		Turn on or turn off buttock air pressure massage, 2 modes options available.
Intensity		2 intensity options available.

## (5) Adjust angle of backrest and footrest

Seat	 Forward	Press this button continually and the seat frame will move forward, loose your finger and the movement will stop.
	 Backward	Press this button continually and the seat frame will move backward, loose your finger and the movement will stop.
Back	 Back Recline	Press this button continually and the backrest frame will move upward, loose your finger and the movement will stop.
	 Back Up	Press this button continually and the backrest frame will move downward, loose your finger and the movement will stop.


## 2.Components install direction

- Place the bottom frame of massage chair on the spa tub.
- Make sure the holes of the bottom frame are aligned with the pre-drilled holes on the spa tub.
- Secure the massage chair onto spa tub by using provided hardware as shown in the diagram below.

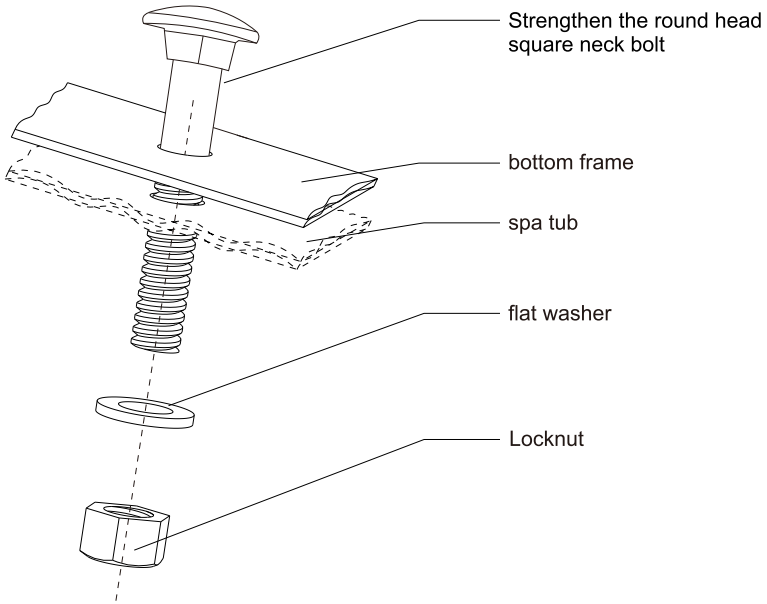


### 3. Armrest assembling instructions

- Insert the position hole of armrest connecting plate into the position pin of armrest support bracket on movable frame.
- Insert hexagon bolt M10x45; elastic washer 10, flat washer 10 sequentially into and threaded through the hole on the armrest supporting seat, screw into the threaded hole on armrest connecting plate, and tighten up.

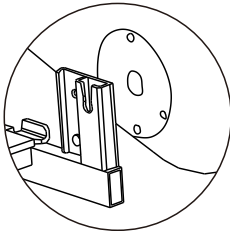


(1)

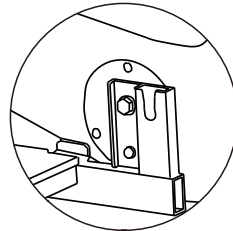


(2)

### 4. Instruction of backrest installation

- Put backrest support shaft into the card slot of backrest support seat. (Note: The shaft can only put into the card slot at certain angle).
- Erect the backrest, the self locking hook will lock the cylinder limit shaft automatically.

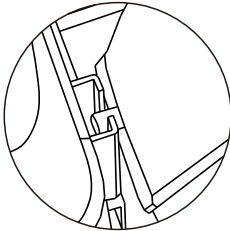

(1)

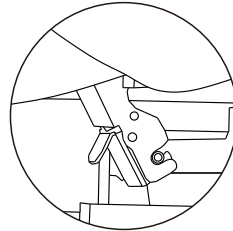


(2)

**Warning:** Before operating the massage, please ensure that the backrest electric cylinder and the pin which fixed backrest connect correctly, and the backrest must be erected before turning on the massage chair. Otherwise, after turned on, the front of electronic cylinder will hit the board under the seat frame, resulting in the seat frame deformation, and even derail or damage the steel ball railway.

### 5. Operation sequence


#### SWITCH ON POWER

- (1).insert the plug into the socket outlet.
- (2).switch on the power button behind the armrest.

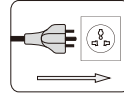


#### MASSAGE

- (1)Select relative massage according to the operation instructions of the control.

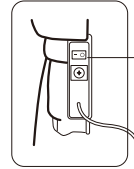

#### SWITCH OFF POWER


- (1).stop overall massage.
- (2).switch off the power button behind the armrest.
- (3).pull the plug out of the socket outlet.

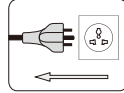

### 6.Attaching or detaching the pillow pad and back cushion



The pillow pad and back cushion (recommended) reduces the kneading intensity of your neck and shoulder, but you can decide to put it on or not according to your personal need.

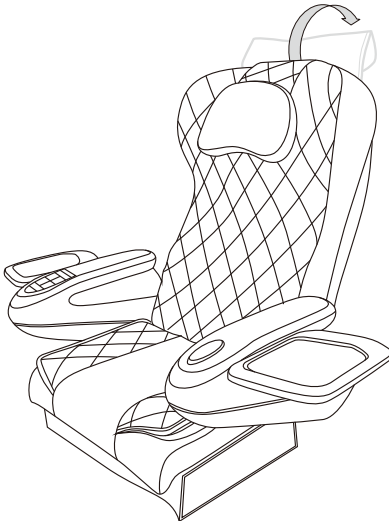**WARNING : IN CASE OF DANGER, DO NOT REMOVE THE CUSHION**




### 7.To Adjust The Back Simultaneously

Backrest adjustment

press and hold the  button on the controller,the backrest will recline backward slowly. Release the button,then the backrest will be stopped upright  
press and hold the  button on the controller,the backrest will move slowly.  
Loose the button,then the backrest will be stopped .

**WARNING: DO NOT EXERT ANY PRESSURE OTHER THAN YOUR OWN BODY ON THE BACK WHILE IT IS RISING.**

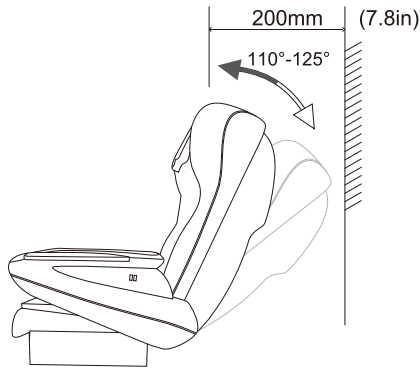


### 8.To Adjust Seat Frame Separately

The flex function of seat frame could fix the different height of the people to adjust the massage place of the calf airbags,the maximum distance moving forward and backward is 120mm.



## Method of Usage

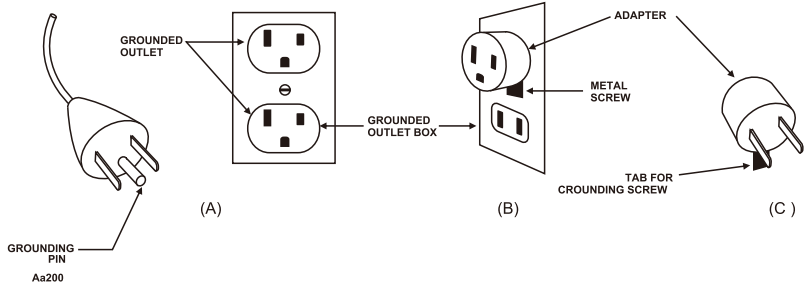



## ■ Product Specification



Model: 9660

Name: Massage Chair Assembly

Rated Voltage: 120V~ 60Hz

Rated Input Power: 100W



Manufacturer of Luxurious Pedicure Spas,  
High Quality Salon and Spa Furniture

**Gulfstream Inc.**

Cambridge, Ontario, Canada N1R 5X5

Tel.: 519-622-0950

Fax: 519-622-6926

[www.gulfstreaminc.com](http://www.gulfstreaminc.com)

**9660 Massage Chair®**, **Gulfstream Inc.®** are trademarks of **Gulfstream Inc.**

Patent #D709619, #8726430, #8680699, #D701318, #8662848, #8531048, #8366418, #D670938, #D670354, #D662340, #D660444, #D650084, #D646082, #D636090, #D632798, #D628704, #D615326, #D613099, #D612509, #D610697, #D600950, #D596883, #D596882, #D590188, #D589285, #D588844, #D583960, #D583576, #D583171, #D565739, #D556480, #7111334, #6991182

© 2021 All rights reserved. Design and Specifications are subject to change without notice.